



Perifusion Chamber Model No. PERI-CHAMBER

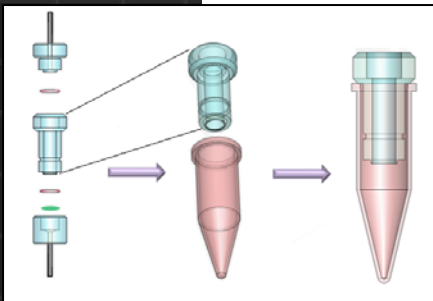

Perifusion Chamber setup

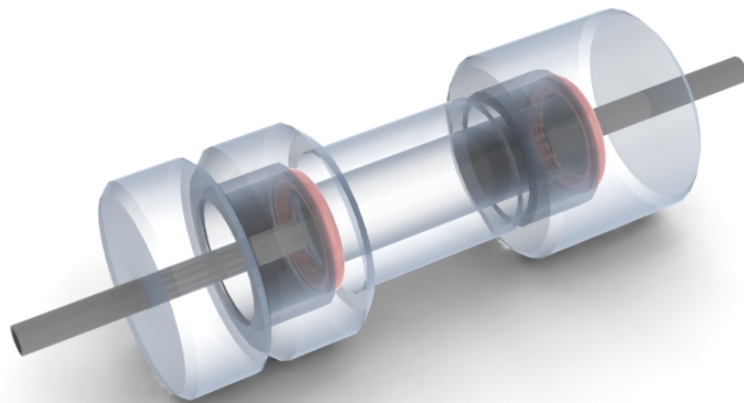

Easy cell recovery

The Perifusion Chamber is specifically designed to use with the Biorep Perifusion System.

It is where all cell reactions will take place. Substances will be delivered through the inlet tube, cells contained in the chamber will react to the substances, and secretions will come out the outlet tube for collection.

- Internal volume = 275 μL
- Optically clear body, polycarbonate material
- Easy cell recovery using a micro-centrifuge tube

*Custom sizes and materials available upon request